

# AutoISF2.2 – Quick User Guide

**Caveat:** *I am not a medically trained person and developed this method purely based on trial and numerical experimentation. For example, I had no mathematical model of the reaction kinetics or of the free fatty acids which are one reason for temporary insulin resistance.*

## Introduction

AutoISF is meant for the advanced user who has a deep understanding of AAPS and who has tuned the system to achieve a TIR of about 90% or better. If such a user is ambitious and wants to improve further, the methods contained can very well help.

This document describes how and when ISF is adapted automatically and provides short descriptions of the weighting factors to tune it to your individual needs. Detailed guides covering special situations and examples will follow later as separate documents including results of using it in full loop mode, i.e. without carb or insulin entry by the user.

The adaptation of ISF is based on special glucose behaviour and results in an adaptation factor just like Autosense. However, here only ISF is adapted and no other setting. The scenarios analysed typically cover the last 10-30 minutes and therefore autoISF reacts much faster to problems or recoveries. Often Autosense would drive me into hypos because of its delayed reaction even after things had come back to target and I disabled Autosense although both can coexist.

AutoISF is part of oref1 in OpenAPS SMB and cannot coexist with the recent DynamicISF which therefore is included in its own plug-in as an alternative to OpenAPS AMA and OpenAPS SMB.

There are 4 different effects in glucose behaviour that autoISF checks and reacts to:

1. **acce\_ISF** is a factor derived from acceleration of glucose levels
2. **bg\_ISF** is a factor derived from the deviation of glucose from target
3. **delta\_ISF** and **pp\_ISF** are factors derived from glucose delta
4. **dura\_ISF** is a factor derived from glucose being stuck at high levels

Finally these factors are compared among each other and against Autosense. Normally the strongest of them will be used with some exceptions as detailed further below. These factors work the same way as the Autosense factor works, i.e. the profile sensitivity is divided by the factor to deliver a final sensitivity ISF.

In the SMB-tab, section Script debug, you can always see which values were assigned during the last loop execution. It also lists explanations in case the factor had to be modified or why it cannot be used. Some interim values like dura\_ISF\_average are listed, too. All that can be seen in the SMB-tab, sections Glucose status, Profile and Script debug, respectively.

Again in analogy to Autosense there are upper (autoISF\_max) and lower (autoISF\_min) limits for how far ISF can be modified in total.

In analogy to enabling SMB there is a setting enable\_autoISF which determines whether any part of autoISF is enabled or none at all.

The settings specific to autoISF are collected in its own menu found at the end of the OpenAPS SMB menu. A screenshot is shown as attachment. Another trick for finding them is to use the filter method at the top of the Preferences page which searches for all settings containing the string entered in the filter field.


## acce\_ISF determination and its impact

This is the latest contribution to autoISF and has been in action since late December 2021.

For things like glucose or delta to change, there first must be an acceleration. Therefore acceleration recognises such changes earlier and is used by autoISF to take pre-emptive action.

Glucose and delta, its 1<sup>st</sup> derivative, play a significant role in AAPS in determining the insulin required. Acceleration, its 2<sup>nd</sup> derivative, was not included so far. One reason might be that it is harder to extract from the glucose history considering that delta needed to be averaged already to provide a reliable signal. In autoISF a best fit algorithm is used to determine the parabola which best matches the glucose data.

Once the formula for the parabola is known it is then very easy to determine the acceleration. Sometimes the fit has bad correlation, i.e. it deviates too much from the glucose readings. In such cases there is no contribution from acceleration and acce\_ISF = 1.

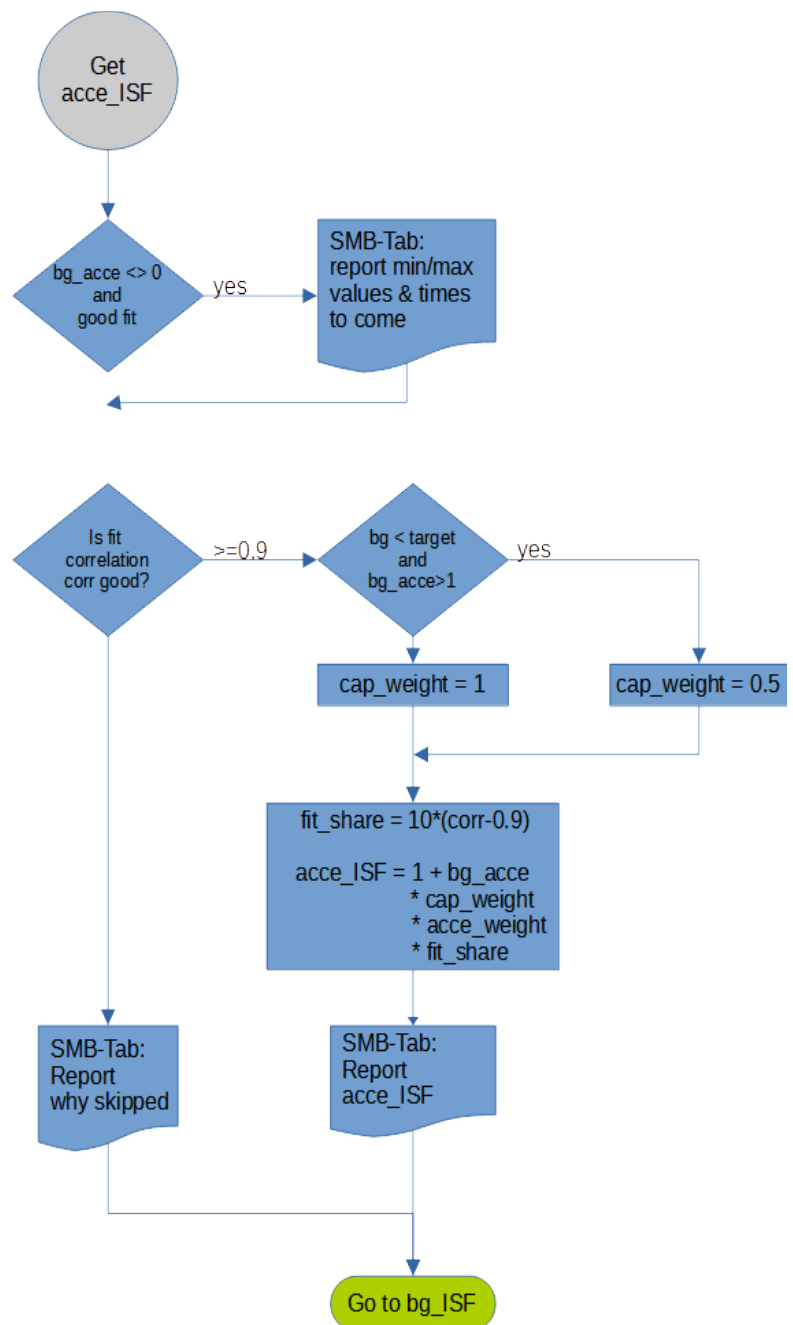

Otherwise acce\_ISF is calculated by


$$\text{acce\_ISF} = 1 + \text{acce\_weight} * \text{fit\_share} * \text{cap\_weight} * \text{acceleration}$$

where   
 fit\_share      a measure of fit quality, is 0% if unacceptable up to 100% if perfect;   
 cap\_weight    is 0.5 below target and 1.0 otherwise;   
 acce\_weight   is bgAccel\_ISF\_weight for positive acceleration   
                  or bgBrake\_ISF\_weight for negative acceleration.


It is still early days but it can be assumed that the weights for accelerating and braking are of similar size. Quite often the acce\_ISF contribution plays the dominant role inside autoISF and is therefore very important and delicate. Weights for acce\_ISF of 0 disable this contribution. Start small with weights like 0.02 and observe the results before increasing them. Keep in mind that negative acceleration will start to happen while glucose is apparently still rising but the slope reduces. Here, acce\_ISF will be <1, i.e. sensitivity grows and less insulin than normal will be required even before the glucose peak is reached.

## bg\_ISF determination and its impact

There are indicators that higher glucose needs stronger ISF. This was evident from all the successful AAPS users defining automation rules which strengthen the profile at higher glucose levels. The drawback is that there are sudden jumps in ISF at switch points and no further or minor adaptations in between.

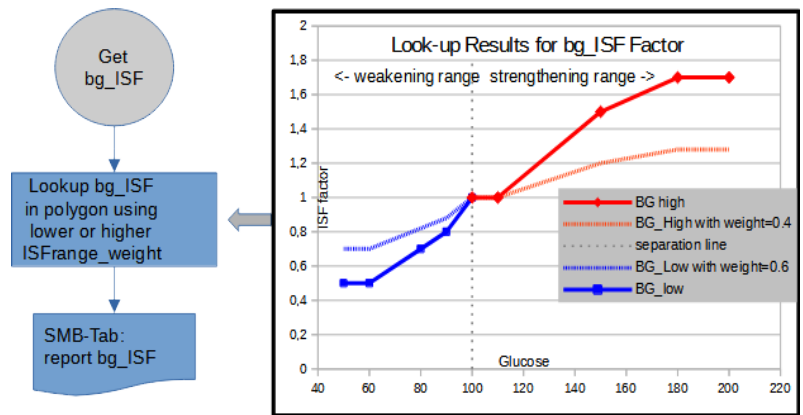

In autoISF a polygon is provided that defines a relationship between glucose and ISF and interpolates in between. This is currently hard coded but the user can apply weights to easily strengthen or weaken it in order to fit personal needs. In principle the polygon itself can be edited and the apk rebuilt if a different shape is required. Developing a GUI for that purpose was considered very tedious especially before knowing whether the results warrant the effort. With this approach you could even approximate the formula well enough that is used in DynamicISF for the ISF dependency.

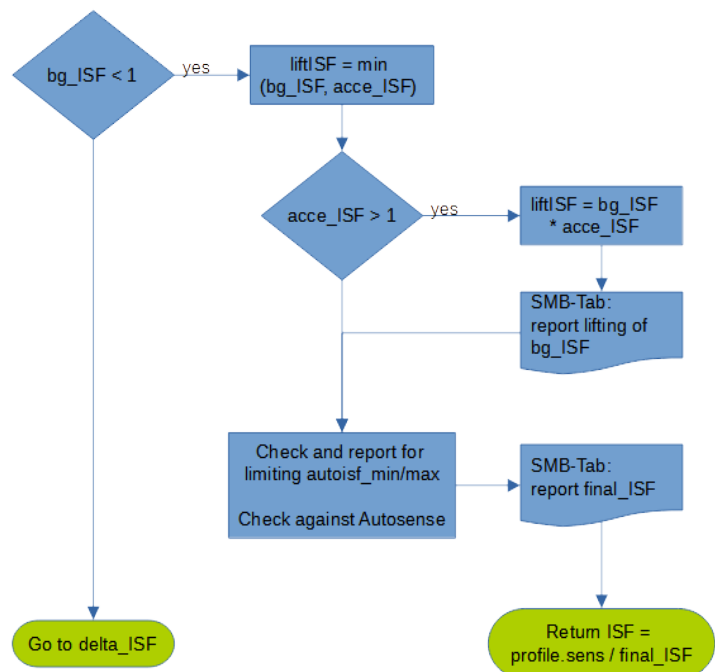

There are two weighting factors depending on whether glucose is below or above target:

*lower\_ISFrage\_weight*

Used below target, weakens ISF the more the higher this weight is; 0 disables this contribution, i.e. ISF is constant in the whole range below target.

This weight is less critical as the loop is probably running at TBR=0 anyway and you can start around 1.

*higher\_ISFrage\_weight*

Used above target, strengthens ISF the more the higher this weight is 0 disables this contribution, i.e. ISF is constant in the whole range above target.

Start with a weight of 0.2 and observe the reactions and check the SMB-tab before you increase it with care.

The result is:

$$bg\_ISF = 1 + xxx\_ISFrage\_weight * glucose\_polygon\_Lookup$$


There is a special case possible, namely below target i.e. when  $bg\_ISF < 1$ . ISF will be weakened and there is no point in checking the remaining effects. Only with positive acceleration the weakening will be less pronounced as that is a sign of rising glucose to come soon.

## pp\_ISF and delta\_ISF determination and its impact

There are two alternative effects in autoISF based on delta. The first one becomes active during meal absorption and was introduced to help users with gastroparesis. They can define a time window (*pp\_ISF\_hours*) during which this effect is active. Alternatively any user can set *enable\_pp\_ISF\_always=true* to activate it all the time. The latter is also the recommended setting for UAM users because in their case no meal start can be detected. Given a positive *short\_avgdelta* and average glucose being above *target+10* the result is:

$$pp\_ISF = 1 + \text{delta} * pp\_ISF\_weight.$$

As a starting value for *pp\_ISF\_weight* use 0.005. Observe the reactions and check the SMB-tab before you increase it with care. A weight of 0 disables this contribution.

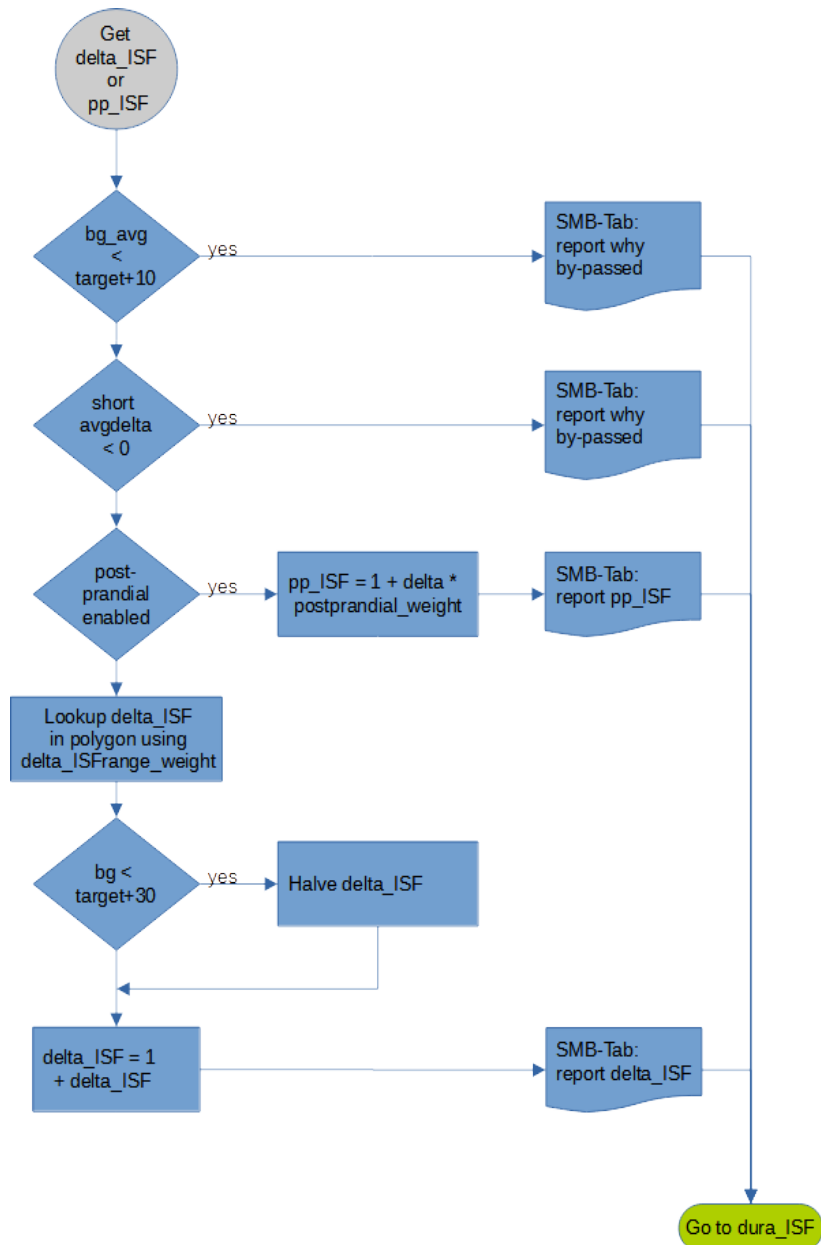

If postprandial settings are not enabled then the alternative method mimics what AAPS users did in defining automation rules that modify profiles or targets based on the level if delta. As in the case of *bg\_ISF* again a polygon is used to lookup an ISF modification. For safety reasons this modification is halved (*target\_penalty*) if *bg<target+30*. Finally the *delta\_ISFrangle\_weight* is applied:

| POLYGON<br>x_data<br>(delta) | POLYGON<br>y_data<br>(Lookup) | delta_ISF<br>with<br>weight=1.0 | delta_ISF<br>with<br>weight=0.06 |
|------------------------------|-------------------------------|---------------------------------|----------------------------------|
| 2                            | 0                             | 1                               | 1,000                            |
| 7                            | 0                             | 1                               | 1,000                            |
| 12                           | 0,4                           | 1,4                             | 1,024                            |
| 16                           | 0,7                           | 1,7                             | 1,042                            |
| 20                           | 0,7                           | 1,7                             | 1,042                            |

$$\text{delta\_ISF} = 1 + \text{delta\_ISFrangle\_weight} * \text{target\_penalty} * \text{delta\_polygon\_Lookup}$$

As a starting value for *delta\_ISFrangle\_weight* use 0.02. Observe the reactions and check the SMB-tab before you increase it with care. A weight of 0 disables this contribution.

## dura\_ISF determination and its impact

This is the original effect of autoISF in action since August 2020. Because autoISF is now a toolbox of several effects this original effect was renamed *dura\_...* It addresses situations when

- glucose is varying within a +/- 5% interval only;
- the average glucose (*dura\_ISF\_average*) within that interval is above target;
- this situation lasted at least for the last 10 minutes (*dura\_ISF\_minutes*).

This is a classical insulin resistance and is typically caused by free fatty acids which grab available insulin before glucose can. Quite often user get impatient in such a situation and administer one or even more rogue boluses. Again and again that leads to hypos later which the *dura\_ISF* approach avoids if carefully tuned.

The method is active if

- no COB is detected like long-time after a meal or in pure UAM mode;
- or *enable\_dura\_ISF\_with\_COB*=true, i.e. even while a meal is still being absorbed.

The strengthening if ISF is stronger the longer the situation lasts and the higher the average glucose is above target:


$$\text{dura\_ISF} = 1 + \frac{\text{avg05} - \text{target\_bg}}{\text{target\_bg}} * \frac{\text{dura05}}{60} * \text{dura\_ISF\_weight}$$

where

$$\begin{aligned} \text{avg05} &= \text{dura\_ISF\_average} \\ \text{dura05} &= \text{dura\_ISF\_minutes} \end{aligned}$$

The user can apply his personal weighting by using *dura\_ISF\_weight*. Start cautiously with a value of 0.2 and be very careful when you approach 1.5 or even higher. By using 0 this effect is disabled.

**CAUTION** for users with Libre CGMs: There can be error states in your CGM which look like long lasting resistance in this context. Therefore Libre CGM users should use *dura\_ISF\_weight* = 0. This way they can still benefit from the other three effects described before.

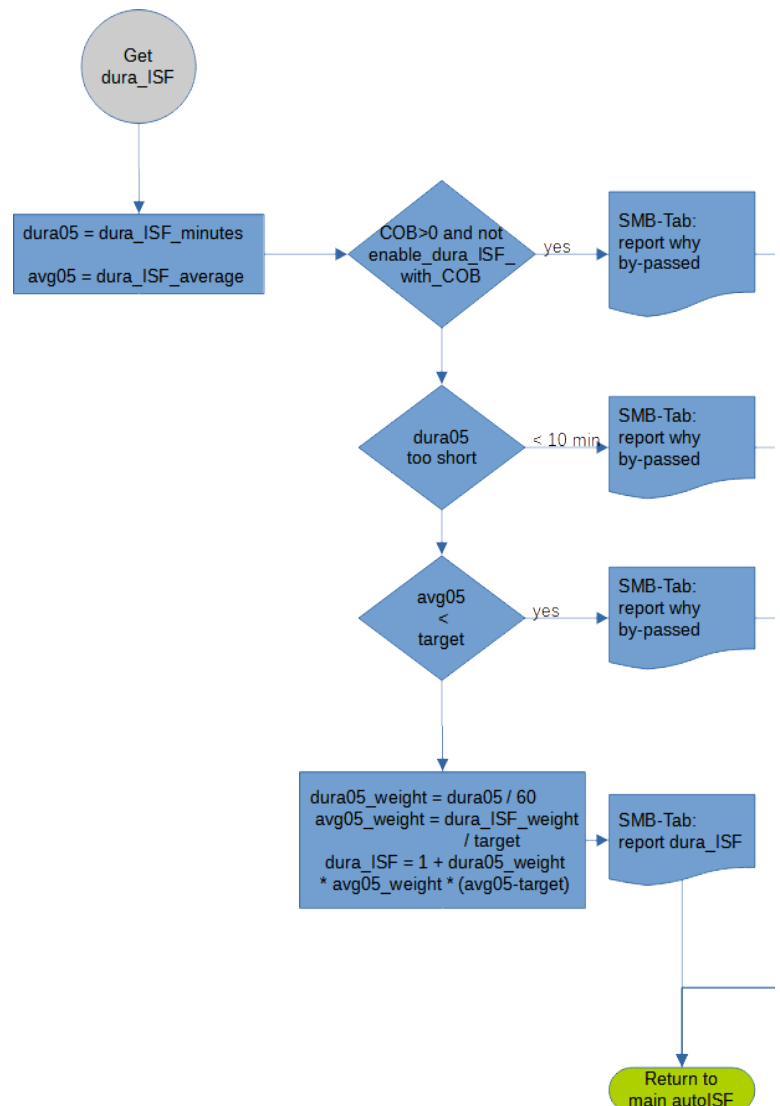


## The combined result from the above factors

Now that the factors for all 4 effects are known how to deduce an end result? The normal case is to pick the strongest factor as the one and only factor to be applied. Here autosense is also part of the game. But how about the exceptions, i.e. when different factors pull in different directions? In order of precedence they are:

- $bg\_ISF < 1$ , i.e. glucose is below target  
If  $acce\_ISF > 1$ , i.e. glucose is accelerating although below target, multiply both factors as a trade-off between them. Then use the weaker of  $bg\_ISF$  and Autosens as the final sensitivity ISF.
- $acce\_ISF < 1$ , i.e. glucose is decelerating while all other effects want to strengthen ISF. In this case the strongest of the remaining, positive factors will be multiplied by  $acce\_ISF$  to reach a compromise. This overall factor will be compared with autosense and the stronger of the two will be used in calculating the final sensitivity ISF.

In all of the above the autoISF limits for maximum and minimum changes will also be applied.

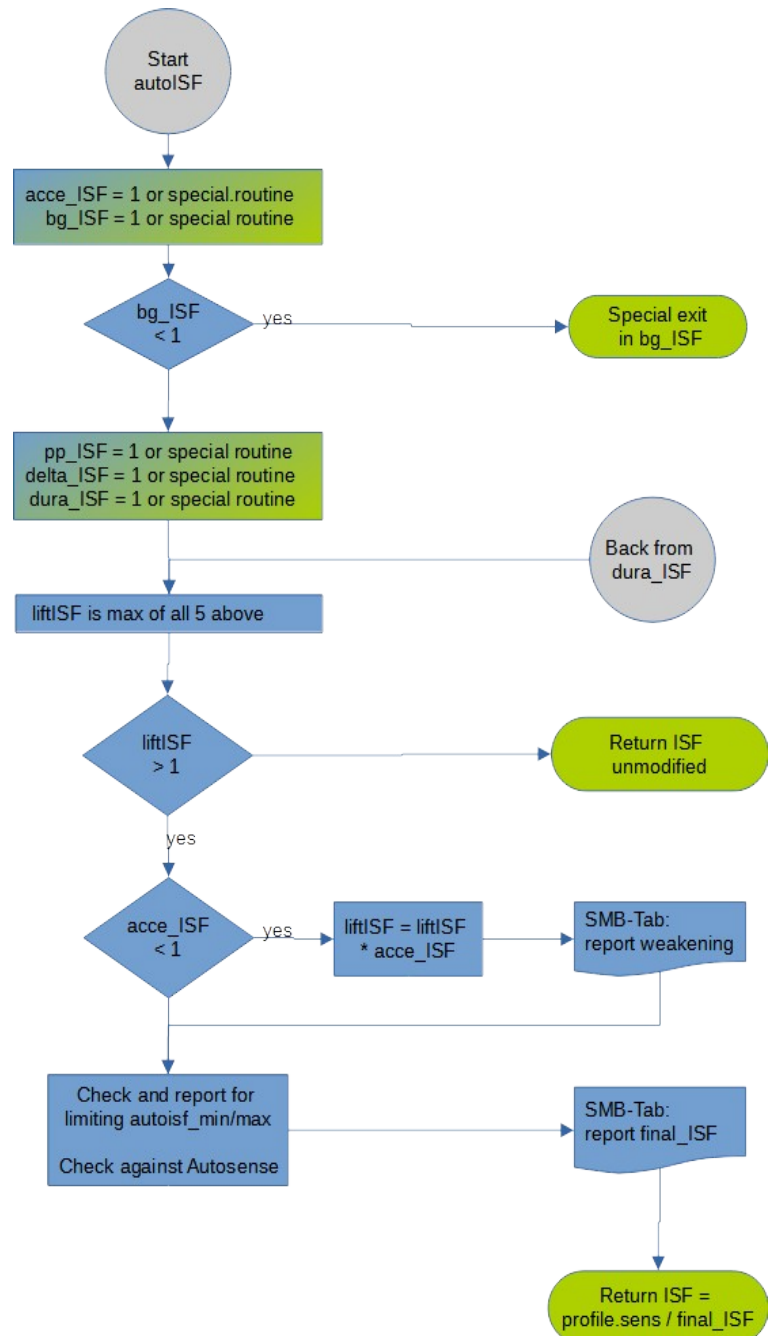



## Adapting SMB delivery

After reading all the methods to strengthen ISF you may wonder whether your maxBolus will always allow as much SMB as requested by the loop with such stronger ISF. This is especially important for pure UAM mode, where you want a few but strong SMBs as soon as meal absorption is detected by acce\_ISF in order to catch up with pre-bolus and meal bolus in standard use of AAPS. Several extensions in autoISF can be used to get there:

- *smb\_delivery\_ratio* is normally hard coded as 0.5 of the insulin requested. This is a safety feature for master/follower setups in case both phones trigger an SMB in the same situation. If this does not apply in your case you may increase this setting to a value above 0.5 and up to even 1.0 if you are very courageous. Best is to leave some margin like the max delivery setting in the calculator.
- Alternatively to a higher but fixed ratio you can use a linearly rising ratio, starting with *smb\_delivery\_ratio\_min* at *target\_bg* and rising to *smb\_delivery\_ratio\_max* at *target\_bg+smb\_delivery\_ratio\_bg\_range*. If *smb\_delivery\_ratio\_bg\_range=0* then this linear rise is disabled and the above *smb\_delivery\_ratio* is used instead.
- *smb\_max\_range\_extension* is a factor that increases the current *maxSMBBasalMinutes* and *maxUAMSMBBasalMinutes* beyond the 120 minute limit.

```
fixed SMB delivery ratio
0.7

variable SMB delivery ratio, lower end
0.65

variable SMB delivery ratio, upper end
0.95

variable SMB delivery ratio, mapped glucose
range
40

SMB/UAM max range extension
3
```

20:17 58%

HOME AKT COMBO **SMB** BEH

UAM Impact: 6 mg/dL per 5m; UAM Duration: 1.2 hours  
minPredBG: 105 minIOBPredBG: 105 minZTGuardBG: 86  
minCOBPredBG: 74  
minUAMPredBG: 108  
avgPredBG: 106 COB: 0.6774242424242392 / 14  
BG projected to remain above 94 for 35 minutes  
BG projected to remain above 77 for 115 minutes  
naive\_eventualBG: 73 bgUndershoot: 4 zeroTempDuration: 115  
zeroTempEffect: 52 carbsReq: -7  
IOB 0.482 > COB 0.6774242424242392; mealInsulinReq = 0.065  
profile mealUAMSMBBasalMinutes: 120 profile current\_basal: 0.4  
**SMB delivery ratio set to interpolated value 0.74**  
naive\_eventualBG 73, sum 0.136/n temp needed, last bolus 0.4m ago,  
maxBolus: 2.4

Result : temp: "absolute"  
bg: 106  
tick: "+3"

## Attachment: Table of all autoISF settings

| Name                        | Use   | Min – Max    | Default | Useful initial value   |
|-----------------------------|---|--------------|---------|------------------------|
| enable_autoISF              | Use any of the autoISF methods  | False - True | False   |                        |
| autoISF_min                 | Lowest ISF factor allowed   | 0.3 – 1.0    | 1       | 0.7 like autosense_min |
| autoISF_max                 | Highest ISF factor allowed  | 1.0 – 3.0    | 1       | 1.2 like autosense_max |
| bgAccel_ISF_weight          | Strength of acce_ISF contribution   | 0.0 – 1.0    | 0       | 0.02                   |
| bgBrake_ISF_weight          | Strength of acce_ISF contribution   | 0.0 – 1.0    | 0       | 0.02                   |
| enable_pp_ISF_always        | Flag to use it even without COB   | False - True | False   | True for UAM mode      |
| pp_ISF_hours                | How long after start of a meal the effect is active                             | 1 – 10       | 3       | 6 for gastroparesis    |
| pp_ISF_weight               | Strength of effect  | 0.0 – 1.0    | 0       | 0.005                  |
| delta_ISFrange_weight       | Strength of effect outside the above postprandial time window                   | 0.0 – 1.0    | 0       | 0.02                   |
| lower_ISFrange_weight       | Strength of bg_ISF effect when bg<target  | 0.0 – 2.0    | 0       | 1                      |
| higher_ISFrange_weight      | Strength of bg_ISF effect when bg>target  | 0.0 – 2.0    | 0       | 0.2                    |
| enable_dura_ISF_with_COB    | Enable dura_ISF even if COB is present  | False - True | False   |                        |
| dura_ISF_weight             | Strength of dura_ISF effect; will also grow with time above target              | 0.0 – 3.0    | 0       | 0.2                    |
| smb_delivery_ratio          | Increase the AAPS standard 0.5 of Insulin required                              | 0.5 – 1.0    | 0.5     | 0.6                    |
| smb_delivery_ratio_min      | Minimum for linearly rising ratio at bg=target_bg                               | 0.5 – 1.0    | 0.5     | 0.5                    |
| smb_delivery_ratio_max      | Maximum for linearly rising ratio at bg=target_bg + smb_delivery_ratio_bg_range | 0.5 – 1.0    | 0.9     | 0.8                    |
| smb_delivery_ratio_bg_range | Width of bg range to reach maximum ratio  | 0.0 – 100    | 0       | 40                     |
| smb_max_range_extension     | Multiplier for maxSMBBasalMinutes and maxUAMSMBBasalMinutes                     | 1.0 – 5.0    | 1       | 1.5                    |



## Attachment: Screenshots of autoISF menu

